

Comfort is king


Corpus®, our most ergonomic and capable power wheelchair seating system, was designed to configure to your body.

Biomechanically optimized power positioning enables functional movements while supporting your optimal position throughout the day.

Dual density backrest foam promotes greater surface area contact while the contoured Ergo seat pelvic well insert can help pressure redistribution and reduce peak pressures.

Stretch-Air® covers offer a breathable, moisture-resistant fabric material to enhance comfort and help protect underlying foam.



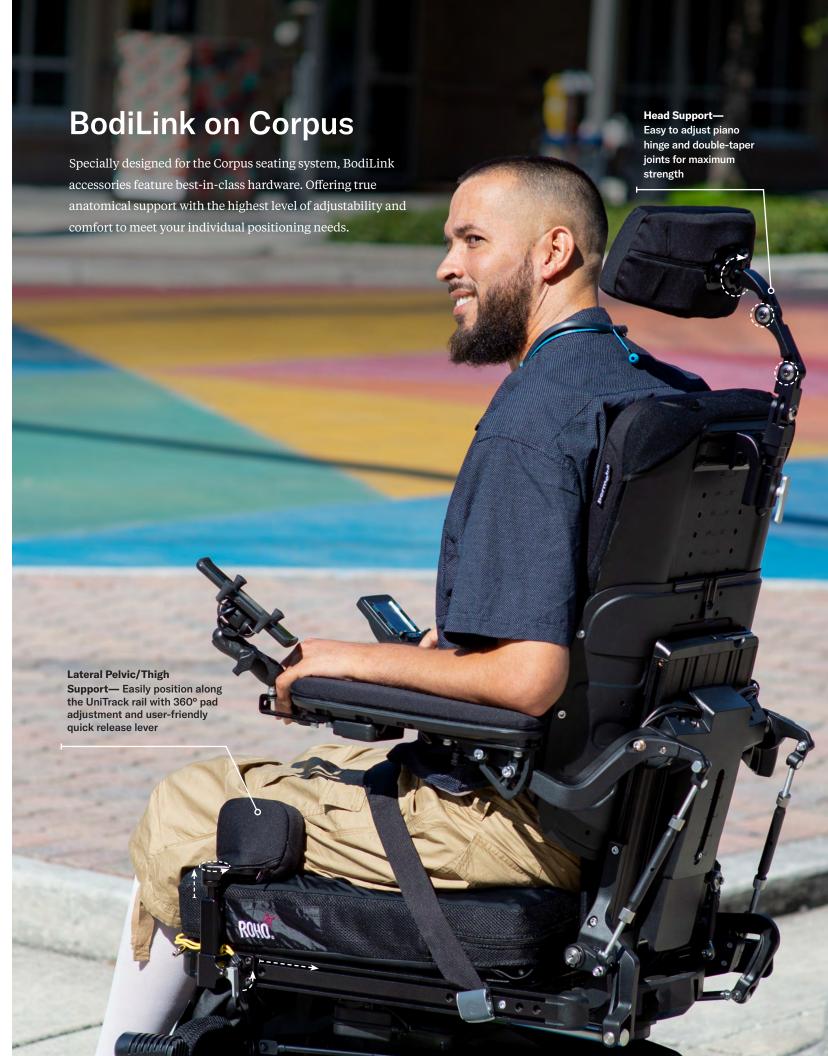





Reach for more

ActiveHeight[™] and ActiveReach[™] combine seamlessly to give you the ultimate advantage in accessing your environment.

Seat elevation with integrated forward tilt creates the most innovative, functional and stable seating solution on the

market today. With up to 14" of vertical seat adjustment and 45° of ActiveReach forward tilt capability, the F-Series empowers you to reach for more and go further faster with new elevated drive speeds of up to 3.2 mph.



Small footprint, big impact

The F3 Corpus offers the smallest footprint of the F-Series front-wheel drive power wheelchairs. Its compact size and agility allows users to easily maneuver in and out of tight spaces. The optimized battery weight distribution, lower center of mass and fully adjustable ComfortRide suspension improve stability while traversing inclines or declines. Experience the enhanced functionality of the F3 by combining outstanding chassis performance with Permobil's industry-leading Corpus Seating System offering 12" of ActiveHeight and up to 30° of ActiveReach forward tilt technology. F3 Corpus Up to 12" ActiveHeight & 30° ActiveReach

36" without support wheels

Strength and stability

Power meets performance on the F5 Corpus with a 7.5 mph top end speed and upgraded ComfortRide Pro suspension. All-wheel independent, fully adjustable, oil dampened shocks provide the ultimate in ride comfort while helping to reduce vibrations that can transfer from terrain up to user. Experience enhanced vibration reduction driving in all seat positions including an unprecedented 14" ActiveHeight seat elevation. Climb and reach higher with increased obstacle negotiation, exceptional motor torque and enhanced stability allowing up to 45° of ActiveReach forward tilt functionality. Standard performance options include aggressive tread tires and caster mud guard protection. F5 Corpus Up to 14" ActiveHeight & 45° ActiveReach Oil-dampened horizontally opposed shocks Caster mud

Stand without compromise

We combined all of the performance features of the F5 Corpus with superior power standing functionality. Fully programmable standing sequences allow user position to be optimized using supine, semi-reclined or sit-to-stand techniques. Further customize standing positions with adjustable chest support, moldable knee blocks and power articulating footplates. From face-to-face interactions to improved health outcomes, the F5 Corpus VS empowers you to stand with confidence.



during standing



Make a statement

With a sleek design, expanded color/hydrographic options and powerful LED lighting, you should be prepared to turn a few heads in your new F-Series.







Hydrographic options

















Liberty

Please note: Actual product colors may vary from colors shown

Platinum

Bronze

Mossy Oak® Breakup[®] Country[™]

Patriot



Specifications

	3 Corpus	F5 Corpus	F5 Corpus VS
Maximum user weight	330 lb	330 lb	300 lb
Maximum speed	6 mph	7.5 mph	5.8 mph, Group 3 7.5 mph, Group 4
Range 13 mi (Gr	o 34) 20 mi (Grp 24) [*]	20 mi [*]	20 mi
Base width	24"	25.5"	25.5"
Race length	th support wheels nout support wheels	42" with support wheels 38" without support wheels	43"
Base turning radius	26.5"	30"	30"
Ground clearance	3"	3"	3"
Recommended battery type G	roup 34 or 24	Group 24	Group 24
Weight including batteries 386 lb (G	rp 34), 407 lb (Grp 24)	412 lb	432 lb
Drive electronics	R-Net 120A	R-Net 120A	R-Net 120A
Seat-to-floor height 17.	5", 18.5", 19.5" [†]	17.5", 18.5", 19.5" [†]	17.5", 18.5", 19.5" [†]
Seat elevation	12"	14"	14"
Drive wheel	14"	14"	14"
Casters	8"	8"	8"
Tilt ontions	terior: 0° – 50° ach: 5°, 10°, 20°, 30°	Posterior: 0° – 50° ActiveReach: 5°, 10°, 20°, 45°	Posterior: 0° – 50° Standing
Legrest angle- power	85° – 180°	85° – 180°	85° – 180°
Backrest angle- power recline	85° – 180°	85° – 180°	85° – 180°
Backrest height 20", 2	23" – 28" (1" incr.)	20", 23" - 28" (1" incr.)	23" – 28" (1" incr.)
Seat depth 14'	- 22" (1" incr.)	14" - 22" (1" incr.)	14" - 22" (1" incr.)
Seat width 17"	- 23" (2" incr.)	17" - 23" (2" incr.)	17" – 23" (2" incr.)

'Actual driving range will vary based on driving and battery conditions. | †Measured without cushion. | Specifications may vary based upon configuration. Please refer to your owners manual for warnings and instructions for the safe operation of your wheelchair in a variety of driving conditions. Chair weights are approximate and includes base, batteries, power seat functions and seating system without accessories.



Note: MyPermobil connectivity is an available option on all new F-Series Corpus & Corpus VS models.





App Store is a trademark of Apple Inc. Google Play is a trademark of Google LLC.

