

Power wheelchairs

F5 Corpus VS




permobil

Stand for more

The F5 Corpus VS that you know and trust is now available with Group 3 and Group 4 base options, making it more accessible than ever to pursue the functional and medical benefits of power standing.

	Group 3	Group 4
Best-in-class base performance	✓	✓
Fully customizable standing sequences	✓	✓
Articulating power footplate	✓	✓
Functional drive speeds in standing	✓	✓
Full range of power tilt, recline, & legrest elevation	✓	✓
14" of ActiveHeight™ seat elevation	✓	✓
Optional no charge accessories: <i>LED light kit, USB charger, MyPermobil</i>	✓	✓
Max speed of 7.5 mph	✗	✓
2-year warranty	✗	✓
Preventative maintenance plan	✗	✓
Hydrographic color options included	✗	✓
One-time shroud color change included	✗	✓

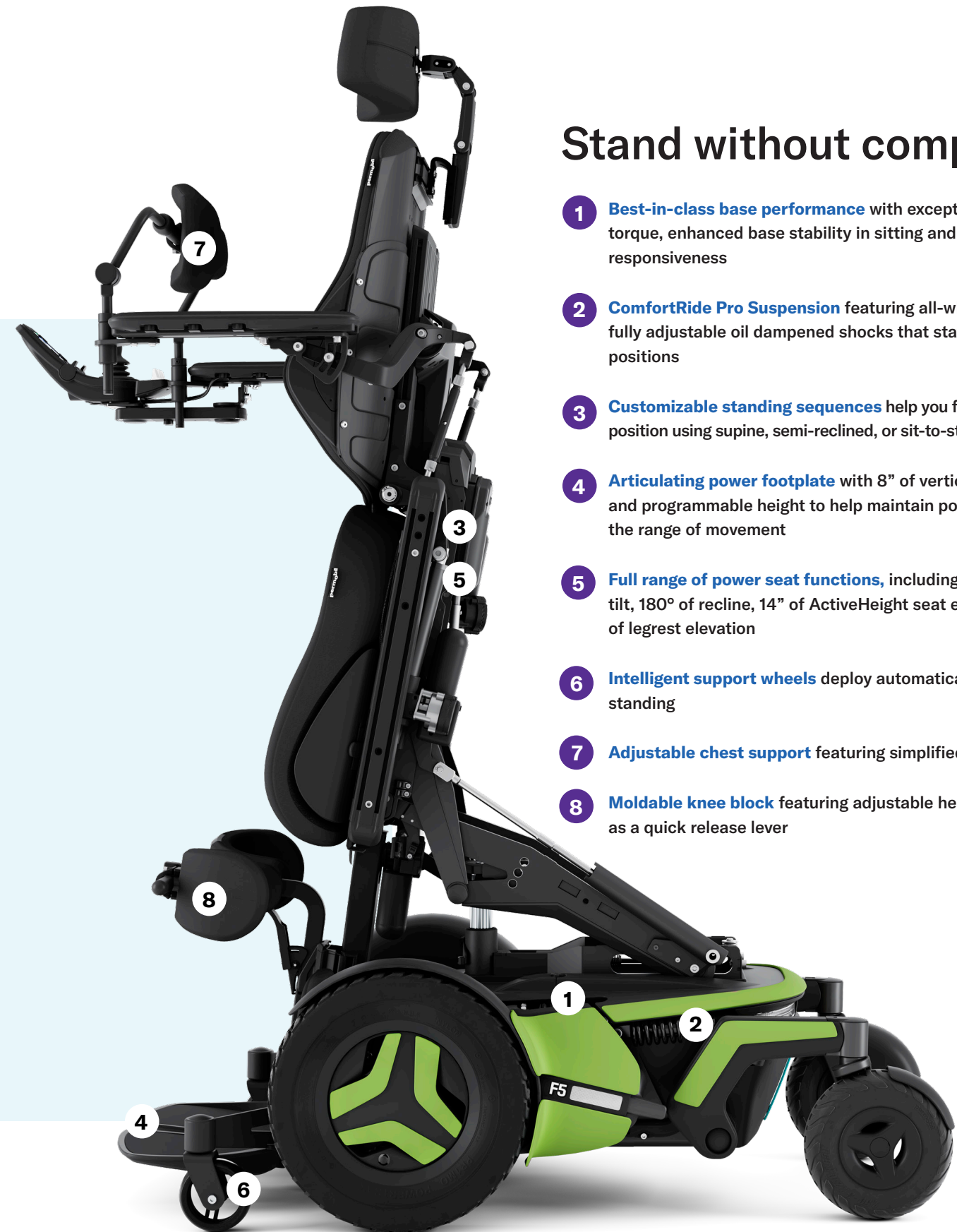

Stand without compromise

- 1 Best-in-class base performance** with exceptional motor torque, enhanced base stability in sitting and standing & chair responsiveness
- 2 ComfortRide Pro Suspension** featuring all-wheel independent fully adjustable oil dampened shocks that stay active in all positions
- 3 Customizable standing sequences** help you find the perfect position using supine, semi-reclined, or sit-to-stand techniques
- 4 Articulating power footplate** with 8" of vertical foot plate travel and programmable height to help maintain positioning throughout the range of movement
- 5 Full range of power seat functions**, including: 50° of posterior tilt, 180° of recline, 14" of ActiveHeight seat elevation, and 180° of legrest elevation
- 6 Intelligent support wheels** deploy automatically for stability in standing
- 7 Adjustable chest support** featuring simplified adjustment points
- 8 Moldable knee block** featuring adjustable height & depth as well as a quick release lever

Your outcome is our top priority

In addition to the superior power standing provided by the F5 Corpus VS, Permobil is proud to offer industry-leading clinical and funding resources. We want to ensure you and your healthcare team have what it takes to achieve an optimal outcome.



Scan the QR code to learn more



Specifications

Maximum user weight	300 lb	Drive wheel	14"
Maximum speed	5.8 mph, Group 3 7.5 mph, Group 4	Casters	8"
Range	20 mi ¹	Tilt options	Posterior: 0° – 50° Standing
Base width	25.5"	Legrest angle - power	85° – 180°
Base length	43"	Backrest angle - power recline	85° – 180°
Base turning radius	30"	Backrest height	23" – 28" (1" incr.)
Ground clearance	3"	Seat depth	14" – 22" (1" incr.)
Recommended battery type	Group 24	Seat width	17" – 23" (2" incr.)
Weight including batteries	432 lb	<i>*Actual driving range will vary based on driving and battery conditions.</i>	
Drive electronics	R-Net 120A	<i>¹Measured without cushion. Specifications may vary based upon configuration.</i>	
Seat-to-floor height	17.5", 18.5", 19.5" [†]	<i>Please refer to your owners manual for warnings and instructions for the safe operation of your wheelchair in a variety of driving conditions. Chair weights are approximate and includes base, batteries, power seat functions and seating system without accessories.</i>	
Seat elevation	14"		